Medical Grade Rapid Anti-viral/Anti-bacterial Decontamination Box (8) UVC Lamps - 304 Stainless Steel - 220V, 50Hz Instruction Manual

Thank you for your purchase of the Larson Electronics IND-TBL-BX-R2 series Decontamination Box.

IMPORTANT

READ CAREFULLY BEFORE INSTALLING THIS FIXTURE. SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE. THE IND-TBL-BX-R2 MUST BE WIRED IN ACCORDANCE WITH NATIONAL ELECTRICAL CODE AND ALL APPLICABLE LOCAL CODES. PROPER GROUNDING IS REQUIRED FOR SAFETY.

WE STRONGLY ENCOURAGE ONLY A LICENSED ELECTRICIAN INSTALL, OPERATE AND MAINTAIN THIS PRODUCT. IF YOU ARE NOT QUALIFIED OR UNFAMILIAR WITH ANY ASPECT OF THIS INSTRUCTION SHEET, CONSULT AN ELECTRICIAN. THERE ARE NO SERVICEABLE PARTS INSIDE.

WARNING: MAKE SURE POWER IS TURNED **OFF** BEFORE STARTING THE INSTALLATION OR PERFORMING ANY MAINTENANCE.


RISK OF FIRE/ELECTRIC SHOCK – DISCONNECT POWER AT BREAKER BEFORE INSTALLING OR SERVICING.
RISK OF PERSONAL INJURY – FIXTURE MAY BECOME UNSTABLE OR DAMAGED IF NOT INSTALLED PROPERLY.
RISK OF BURN – ALLOW FIXTURE TO COOL BEFORE SERVICING.
WARNING: ULTRAVIOLET LIGHT IS HARMFUL TO EYES AND SKIN. DO NOT STARE DIRECTLY INTO THE ULTRAVIOLET BEAM WITHOUT WEARING PROTECTIVE GEAR TO PROTECT EYES FROM LONG TERM DAMAGE. UNPLUG LAMP BEFORE SERVICING.

The Larson Electronics IND-TBL-BX-R2 series is a Medical Grade Rapid Anti-viral/Anti-bacterial Decontamination Box which can disinfect kill 99% of viruses, bacteria, mold, and spores and bacteria found on tools, medical equipment and work devices. With access to 32.5-watt UV lamps, this unit offers UV-C output at a range of 200nm - 280nm (253.7nm). The 304 stainless steel disinfection box comes with a timer, external indicator light and an automatic shut off safety feature.

Wiring

This UV oven is equipped with 25 feet of 16/3 chemical and abrasion resistant SOOW cord that is fitted with an industrial grade cord cap for easy connection to common wall outlets. Alternatively, we also offer flying leads with no cord cap for hard wire applications, or users who wish to install their own cord cap. See below for flying leads instruction. For options with a cord cap installed, simply mate this cord cap with a matching receptacle outputting voltage within the rated range. Check product nameplate for correct voltage.



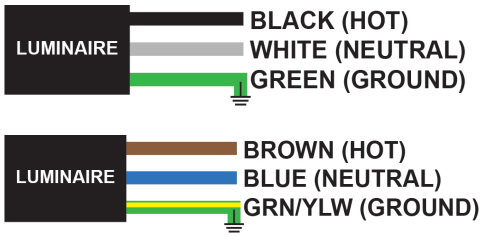


We strongly encourage a licensed electrician install this product, in all locations especially in outdoor areas where weatherproofing may be required.

Warning: *Check product label for correct input voltage!*

Attach supply line wires to the appropriate light fixture wires as shown in the diagram ➔

Secure each pair of wires according to the diagram. Utilize a junction box or similar device and take precautionary steps for weatherproofing all connections if installed in a location where water may come in contact with the unit. Ensure the unit is properly grounded and that wiring is done according to all local and national electrical/building codes.

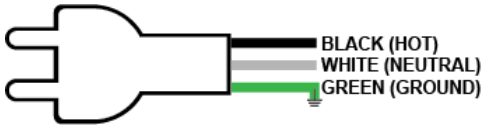

Bulb Installation and Replacement

Warning! - Do not touch the bulbs with your bare hands! The oils in your skin will reduce performance and significantly reduce the lifespan of the bulb unless removed with a soft, non-abrasive cloth and alcohol.

Your decontamination box is shipped with the bulbs removed to prevent damage during shipping. You must install the bulbs before operating. In order to achieve this, you must first remove the protective cage that will cover the bottom bulbs inside the box by removing the (4) screws securing it in place, one at each corner. Once the cage is removed, begin installing the bulbs in any order by lining up the 2-pin connector at the ends of the bulb with the slots in the bulb receptacles, insert, and

turn 1/4 turn to secure the bulb in place. Repeat this process for the remaining bulbs, and then reinstall the protective cage in the bottom of the decontamination box. Follow these same instructions for changing the bulbs, should the need arise. Only use bulbs of the same type and rating when replacing.

Operation

Once your decontamination box is connected to power and bulbs are installed, simply place the unit on a flat surface capable of supporting the load. Unlatch the front facing door, and swing it downward to the open position. Verify that the object(s) to be sanitized will fit inside and clear the bulbs. Place your objects to be sanitized inside, close and latch the door, set the timer to the desired decontamination time, and press the green start button. Use the orange knob on the side of the box to set the timer. An external indicator light notifies the user that the light is on and sanitation is ongoing. Once the timer duration has elapsed, the light will automatically power off. There is also a pressure switch located inside the door frame of the box that will turn the lights off if the door is opened during operation, preventing exposure to harmful UV rays.

USE AND CARE

Unauthorized modification may impair the function and/or safety of this device and could affect the life of the equipment. Always check for damaged or worn out parts before using the device. Store it in a secure place out of the reach of children when not in use. Inspect for good working condition prior to storage and before re-use.

REPLACEMENT PARTS

The IND-TBL-BX-R2 is designed to provide years of reliable performance. Should the need for replacement parts arise, please contact Larson Electronics.

THESE INSTRUCTIONS MAY NOT COVER ALL DETAILS OR VARIATIONS OF THIS PRODUCT FOR YOUR EQUIPMENT OR INSTALLATION REQUIREMENTS. SHOULD FURTHER INFORMATION NOT COVERED BY THESE INSTRUCTIONS BE REQUIRED, PLEASE CONTACT LARSON ELECTRONICS BY EMAIL AT SALES@LARSONELECTRONICS.COM OR BY PHONE AT 1-800-369-6671 FOR FURTHER ASSISTANCE.

PLEASE VISIT LARSONELECTRONICS.COM FOR **WARRANTY** AND **RETURN** INFORMATION.

Ultraviolet Radiation Safety

WARNING: Do not attempt operation until you are familiar with all warnings, precautions, and procedures outlined within this instruction sheet.

OPERATION PRECAUTIONS

- **UV LIGHT WARNING:** UV-C light can cause temporary or permanent loss of vision, and temporary acute redness or ulceration (mild to severe sunburn) to exposed skin. To prevent exposure, do not operate the device in any application that allows UV-C light to be visible during operation.
- **UV LIGHT CAUTION:** UV light may degrade plastic and rubber components after longterm exposure.
- Currently, there are no work place related rules and regulations that are set by OSHA (Occupational Safety and Health Association) in regard to UVC environmental health and safety.
- **NEVER LOOK DIRECTLY AT THE BULB OR DISINFECTING TARGET ZONE**

OPERATION SAFETY

- People and animals are forbidden from being exposed to the light output of this light for extended periods of time. Long term exposure will case damage to the eyes and skin.
- This product produces UV radiation. Anything in front of the unit will be exposed to UV radiation. Operators positioning the light while not illuminated will not be exposed. Radiation is only emitted while the light is turned on. Radiation will not linger in the environment when the lamps have been powered off. Safety equipment is only required when the device is in-use.
- Wear eye protective gear when operating the equipment. Polycarbonate safety glasses at minimum are required during all operations. Full frontal polycarbonate or UV-preventative face masks are recommended.
- Wear long sleeves and gloves when operating the equipment. Personal protective gear is required for long term exposure.

OPERATION SAFETY STRATEGIES

- Instruct service personnel to never look directly at UV-C light without adequate eye protection.
- Make it a policy to never enter areas being disinfected with UV lighting systems without powering the system down first. Methods of powering down depends on the options of the device purchased and may include unplugging the light source, operating the power switch to power the lamps off, or turning the system off via a wireless remote control or wireless application.
- Place warning labels near all openings/entrances to areas being disinfected with UV lighting systems. When possible, lock access doors to prevent unexpected guests from entering the area.
- Close blinds/curtains of areas being disinfected with UV lighting systems if available. If not available, block off direct access to windows/doorways to prevent access to the area. If you cannot fully block off all access (hallways with visible windows, etc) place exterior warnings in areas prior to the exposed windows (hall pedestals, banners, posters) explaining what areas are being disinfected and to not look into the room.
- Educate non service personnel about the intended operations for UV disinfecting and what signage will be posted. Instruct non service personnel to never look directly at UV-C light.
- Instruct service personnel to always wear personal protective gear (as outlined above) when they must operate in an area that has UV-C exposure, for short or long term work.

These safety precautions will help ensure service personnel are protected from accidental exposure while maintaining the effectiveness of UV-C to eradicate biological contaminants.

In the event of UV exposure, the following actions are recommended.

- See an ophthalmologist if eye damage is suspected.
- Treat skin lesions immediately.
- For severe skin lesions, seek local medical treatment and follow their recommendations for treatment.
- Follow your organization's incident reporting procedure. These often require documentation of the date and time of the incident, persons involved, equipment involved and type of injury.

THESE INSTRUCTIONS MAY NOT COVER ALL DETAILS OR VARIATIONS OF THIS PRODUCT FOR YOUR EQUIPMENT OR INSTALLATION REQUIREMENTS. SHOULD FURTHER INFORMATION NOT COVERED BY THESE INSTRUCTIONS BE REQUIRED, PLEASE CONTACT LARSON ELECTRONICS BY EMAIL AT SALES@LARSONELECTRONICS.COM OR BY PHONE AT 1-800-369-6671 FOR FURTHER ASSISTANCE.

PLEASE VISIT LARSONELECTRONICS.COM FOR **WARRANTY** AND **RETURN** INFORMATION.